



Designed to meet the special needs of your patients weighing up to 600 pounds, the Invacare® Reliant™ 600 heavy-duty lift is secure and sturdy, providing a safe alternative to manual lifting. Both durable and functional, the Reliant 600 is also a reassuring and comfortable patient lift. When answering the needs of every patient is your priority, this is your wisest choice.



Model no. RPL600-1 Power lift with low base.

Features

- High lift range allows transfer to high surfaces
- Low friction casters significantly improve rollability and contribute to caregiver safety
- Designed to meet the needs of the growing bariatric population up to 600 lb.
- Manual lowering for instant response to power loss
- Unique sling fabric designed to provide the utmost in comfort and strength
- Available with power opening base (model no. RPL600-2)
- Essential device for safe patient handling in hospitals and LTC (long-term care) facilities
- Limited warranty includes three years on lift, one year on electric components
- Rear locking casters are standard
- Manufactured in an ISO13485 facility

Lift Specifications

	Power Adjustable RPA600-1	Power Low Base RPL600-1	Power Lift with Pwr.Opening Base RPL600-2
Height			
Maximum at sling hook-up	70"	68"	68"
Minimum at sling hook-up	29.5"	28"	28"
Base height (clearance)	6.8"	4.5"	4.5"
Base width			
Open	41"	41"	41"
Closed	26.5"	26.5"	26.5"
Base length	48"	48"	48"
Casters			
Front	5"	3.5"	3.5"
Rear	5"	5"	5"
Weight capacity	600 lb.	600 lb.	600 lb.
Product weight			
Out of carton	112 lb.	112 lb.	115 lb.
Shipping	138 lb.	138 lb.	141 lb.



Power lift with power opening low base.


Model no. RPL600-2



Invacare® Reliant™ Digital Scale

Model no. RLS6

- Precision scale mounted on swivel bar accurately measures weights up to 600 lb.
- Retrofits on all Reliant™ 450/600 lifts. Model 9805 requires Reliant swivel bar
- Easy-to-operate, four-button display
- Automatic shut-off feature prolongs battery life
- Two-year limited warranty

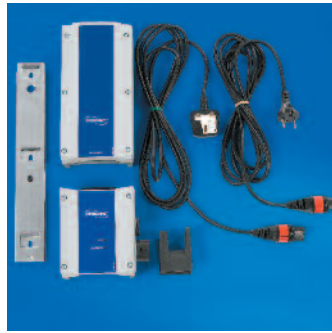


Charging Station

- One 24-volt battery on mast with quick release
- Battery can be replaced immediately with second battery in wall-mounted charger
- UL listed and CSA certified




Manual Lowering

- Manual emergency lowering feature gently and safely lowers resident in case of total power loss
- Essential safety feature
- Manual lowering is easy and effortless. Simply insert the end of a ball point pen into the hole marked "emergency"; lift on the accessible ring-pull

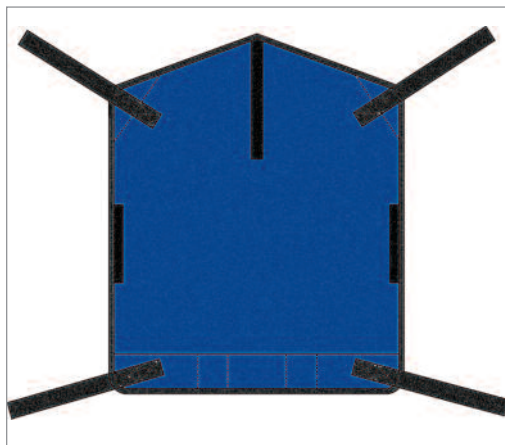

Model no. R140 Mesh sling

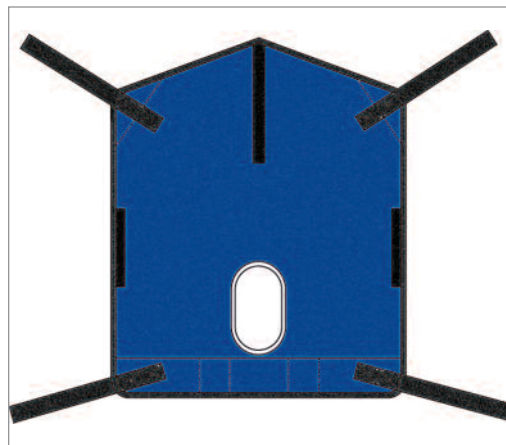

Model no. R141 Mesh sling with commode opening

Sling Specifications

	R140	R141
Width	45"	45"
Commode opening	n/a	7"
Length	55"	55"
Commode opening	n/a	13"
Weight capacity	600 lb.	600 lb.
Product weight		
Out-of-carton	2.2 lb.	2.1 lb.
Shipping	2.7 lb.	2.6 lb.

Complete instructions for lift usage and sling procedures are found in the Invacare Operating Manual. Only Invacare slings and accessories should be used on Invacare lifts.